

MEIKO TOPIC

Washer-disinfectors





MORE TIME FOR PEOPLE

Just a push of a button and a MEIKO washer-disinfector will clean and disinfect your care utensils. Bacteria and viruses are effectively reduced in accordance with DIN EN ISO 15883 – even C. diff spores!

This means: hygiene safety on the one hand, and a reduced workload on the other hand.

Even monitoring can shine! Jörg Menzel, Managing Director of Blindenwohnstätten (housing for the blind) in Berlin agrees: 'In caregiving, we need technology that lightens our workload and makes processes more efficient.'

Depending on what is needed in the space, *TOPIC* is available floor mounted, as a table top unit or with ground clearance to mount on the wall. We can even provide it as a fully equipped combined care unit, containing everything you may need.

TOPIC individual appliances are specially designed to quickly and easily replace out-of-date technology. They are modified to fit into the position of the old appliance, whatever the make of the original – and to use the same connections.



WHAT TO EXPECT:

- Hygienic safety due to A_0 value control system to DIN EN ISO 15883
- Steamtight cleaning and disinfection system
- 3 powerful programme cycles (short, normal, intensive programme)
- Ideal handling with the 'One Touch' operating concept
- Lower costs thanks to reduced use of operating resources
- Flexibility thanks to multifunctional holders
- All-round carefree package: the *MEIKO Clean Solution Circle* including chemicals, service, academy & planning

TOP TECHNOLOGY

The secret to its durability? Solid technology made of stainless steel – for a solid investment.

Top: the deep-drawn wash chamber has no corners or tight angles. It is designed to be ideal for self-cleaning. The ultra-efficiently positioned power nozzle system ensures optimal cleaning of urine bottles, bedpans, commode buckets etc.

You can direct the programme effortlessly using our one-touch technology. From starting a cycle to carrying out an analysis – this is all done using the membrane keypad.

MORE TOP FEATURES:

1. Completely steamtight system


Benefits

- No danger of germ escape
- High workplace and hygiene safety
- Optimised moisture protection

2. A₀ values checked by temperature sensors that monitor each other

Benefits

- High level of hygiene safety (according to EN ISO 15883-1 and -3*)
- Monitoring of disinfection performance and the cleaning process

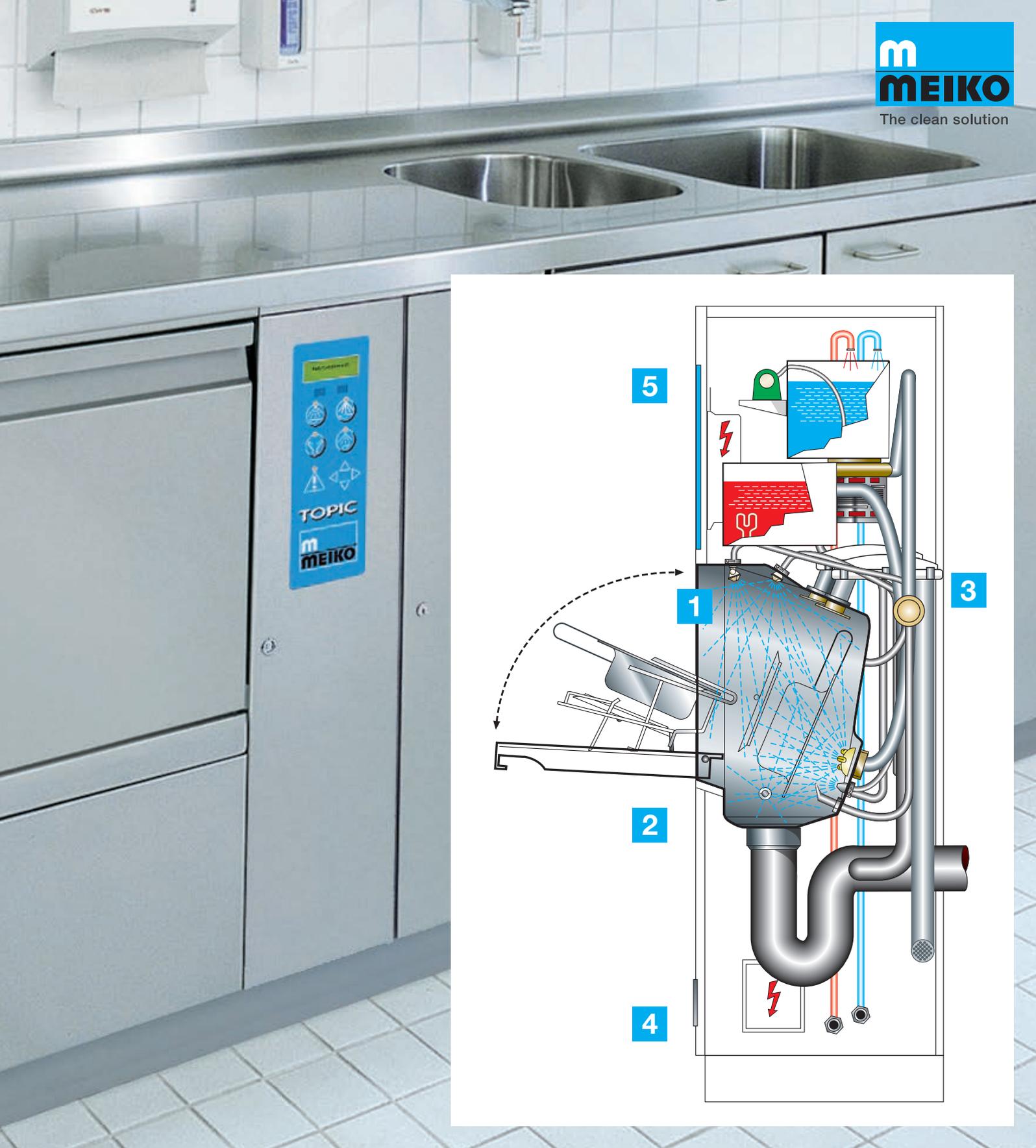

3. System disinfection of all water-carrying lines after each batch cycle

Benefits

- No re-contamination of the system
- Consistent steam distribution and disinfection performance on all water-bearing components

4. Unambiguous coding for detergent (optional) as well as rinse aid and water softener

Benefits

- Maintains process reliability
- User safety
- Makes work easier for staff

5. Sensor-controlled monitoring of dosing of water softener with rinse aid and detergent

Benefits

- Excellent operational safety
- Optimal use of operating resources
- Prevents limescale build-up

MEIKO TOPIC – machine upgrades



RECOOLING WITH DRYING EFFECT (optional)

- For cooling and drying care utensils with filtered air
- Water savings of up to 8 litres per programme cycle (compared to recooling with water)

ADJUSTABLE A₀ VALUE:

- For flexible setting of the A₀ value
- The available settings: A₀ 60, A₀ 600, A₀ 3,000 (3,000 only in connection with recooling and drying)
- Based on the recommendations of your infection control team for basic hygiene and as needed

AUTOMATIC DOOR WITH INFRARED HAND SENSOR AND/OR FOOT SWITCH: (optional)

- Minimises the risk of germ escape
- Simplifies handling in caregiving

THE PERFECT A₀ VALUE IN THREE STEPS



Step 1:

Choose the cleaning programme based on the type of care utensils and the degree of soiling – press the button and hold it down.



Step 2:

Display shows the available A₀-settings: A₀ 60 / A₀ 600 / A₀ 3,000. Release the button when the setting you want appears in the middle.



Step 3:


The programme starts automatically as soon as you release the button.

CLEAN SOLUTION CIRCLE

We will support you at every stage – from planning services through to machine installation and maintenance and beyond to staff induction and training with the necessary chemical products – ensuring that you never have to compromise on hygiene.

Your one-stop provider:

The *MEIKO Clean Solution Circle* brings in all relevant parameters for the cleaning process and their relationships and dependencies on one another. At its core is ‘Sinner’s circle’ of time, temperature, mechanical action and chemicals – the four factors that chemist Herbert Sinner identified as the key to cleaning performance. Yet only with the right consulting and use, the right services and workflows, and training from the *MEIKO Academy* does the complete, holistic picture come into view. Then hygiene, environmental considerations and economy are assured.

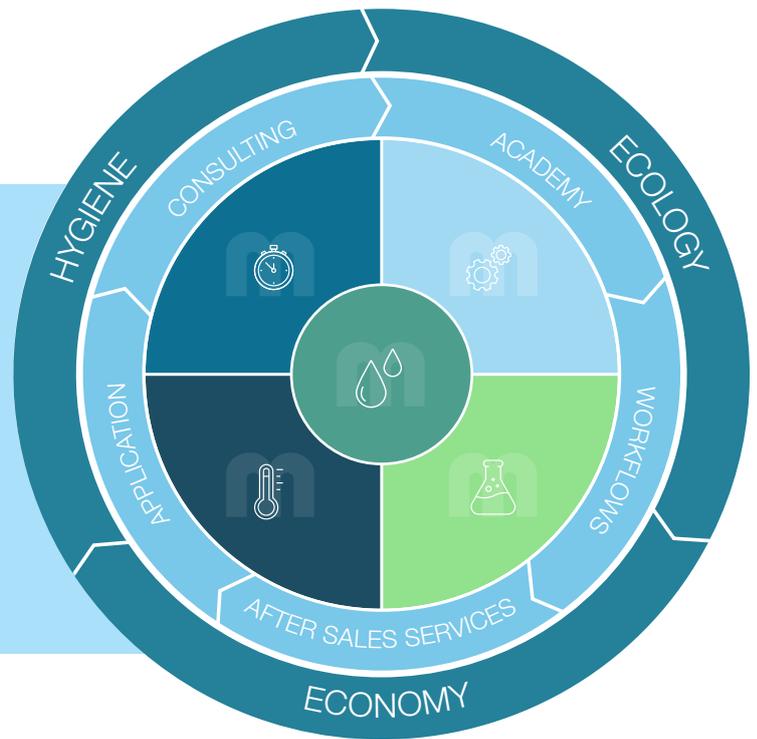

MEIKO CHEMICALS

To achieve perfect cleaning results, the technology, detergent (optional) and rinse aid have to be suited to one another. That is exactly what happens with *DOYEN* products approved by MEIKO.

Benefits of *DOYEN* chemicals:

- Approved for use in *TOPIC* appliances
- Optimum dosing quantities which are kept as small as possible
- Tested for compatibility with materials



MEIKO TOPIC – all-round carefree package



MEIKO PLANNING

Whether an individual appliance or a combined care unit: we plan your complete utility room together with you. That's how we optimise hygienic utility room workspaces. The results meet your needs and fulfil all statutory requirements.

Our brand distinctives include:

- Clean solutions, made to measure. For refits, renovations or new builds
- Seamless connections
- Sealed surfaces to perfectly connect the individual components
- All stainless-steel



MEIKO ACADEMY

Did you know that we also offer training for in-house technicians? The *MEIKO Academy* provides direct access to the knowledge base that we have built up over the decades. We administer online and in-person courses to put you on the cutting edge of technology. Our training courses are rounded off with important information from the world of microbiology, infection prevention and legal aspects.



Find out more about our training courses at:
www.meiko-academy.com

MEIKO SERVICE

Our service means: fast response times, short journeys, high availability of spare parts in the

service vehicle. The advantages of *MEIKO Service*: 365/24/7 service hotline, authorised service technicians (*MEIKO Academy* trained).



MEIKO TOPIC – model overview

Individual appliances

MEIKO TOPIC 10 – wall-mounted unit (with or without slop sink unit)

For off-the-floor wall mounting for easy floor cleaning



MEIKO TOPIC 10

Dimensions (mm): W 500, D 450 or 588, H (appliance) 1,630



MEIKO TOPIC 10 A

Dimensions (mm): W 1,000, D 450 or 588, H (appliance) 1,630
Splash guard (optional)

MEIKO TOPIC 20 – floor-standing unit (with or without slop sink unit)
For floor mounting with appliance base for fast installation



MEIKO TOPIC 20
Dimensions (mm): W 500, D 450 or 588, H 1,730



MEIKO TOPIC 20 A
Dimensions (mm): W 1,000, D 450 or 588, H 1,730
Splash guard (optional)

MEIKO TOPIC – model overview

Individual appliances

MEIKO TOPIC 40 – table top unit (with or without hand wash basin)

For floor mounting with appliance base for expanding the work and storage surface

- Equipped as standard with a lockable storage cabinet for chemical containers
- Tabletop available with upstand for an additional charge



MEIKO TOPIC 40
Dimensions (mm): W 900, D 600, H 900



MEIKO TOPIC 40 B
Dimensions (mm): W 900, D 600, H 1,000
Hand wash basin: B 340, T 370, H 150

Combined care units

MEIKO combined care units are available with floor-standing units or table top units in several versions to suit the conditions in the room and customer requirements. The units provide a multitude of functions that make everyday care work a lot easier.

SAN 14 B / MEIKO TOPIC

floor-standing or table top unit with slop sink unit, hand wash basin

Dimensions (mm): W 1,400, D 600, H 1,000
Slop sink unit: H 700
Slop sink basin: Ø 362
Hand wash basin: W 340, D 370, H 150



SAN 19B – floor-standing or table top unit with slop sink unit, hand wash basin and cabinet (500 mm)

Dimensions (mm): W 1,900, D 600, H 900
Slop sink unit: H 700
Slop sink basin: Ø 362
Hand wash basin: W 340, D 370, H 150
Cabinet space: W 500



SAN 20BW – floor-standing or table top unit with slop sink unit, hand wash basin, wash basin and cabinet (600 mm)

Dimensions (mm): W 2,000, D 600, H 1,000
Slop sink unit: H 700
Slop sink basin: Ø 362
Hand wash basin: W 340, D 370, H 150
Wash basin: W 500, D 400, H 300
Cabinet space: W 600



SAN 24W – floor-standing or table top unit with slop sink unit, wash basin and cabinet (1,000 mm)

Dimensions (mm): W 2,400, D 600, H 900
Slop sink unit: H 700
Slop sink basin: Ø 362
Wash basin: W 500, D 400, H 300
Cabinet space: W 1,000



SAN 29BW – floor-standing or table top unit with slop sink unit, hand wash basin, wash basin and cabinet (1,500 mm)

Dimensions (mm): W 2,900, D 600, H 900
Slop sink unit: H 700
Slop sink basin: Ø 362
Hand wash basin: W 340, D 370, H 150
Wash basin: W 500, D 400, H 300
Cabinet space: W 1,500



MEIKO TOPIC – accessories

In addition to our MEIKO cleaning and disinfection technology, we also offer a wide range of accessories and furnishings. That means we can enable safe and ergonomic work patterns, as well as hygienic storage for your care utensils. We would be happy to develop a range of suggested solutions especially for you!



Drying and storage racks*

- For disinfected urine bottles, bedpans and containers with a storage grid or drip tray
- Available in widths from 500 to 1,400 mm



Wall shelving*

- With height-adjustable shelves
- In all widths (depth 350 mm)



Closed wall-mounted cabinets*

- Cupboard units with hinged doors and one intermediate shelf
- Available in widths 500/600/800/1,000/1,200/1,500 mm
- Optional: light screen and bulb



Open wall-mounted cupboards

- Cupboards for open storage with one intermediate shelf
- Available in widths 500/1,000/1,200 mm



Cleaning combination*

- Base cabinets with drawers and other features for combined care units
- In a variety of widths and heights




Slop sink unit*


- Standard model optionally with built-in slop sink flushing
- Can be controlled separately from the washer-disinfector

* Pictures contain special equipment


** Optional: light screen and bulb

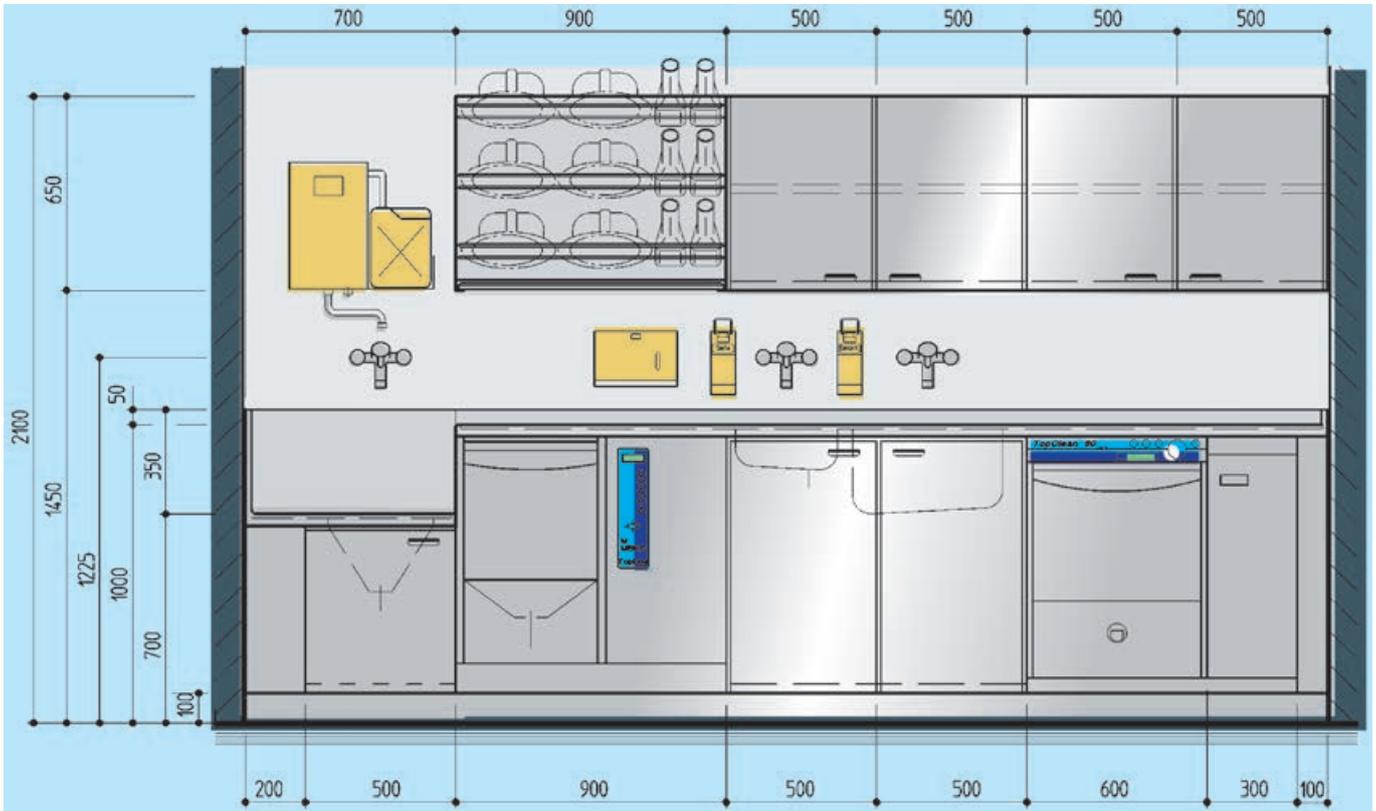


MEIKO TOPIC floor-standing unit

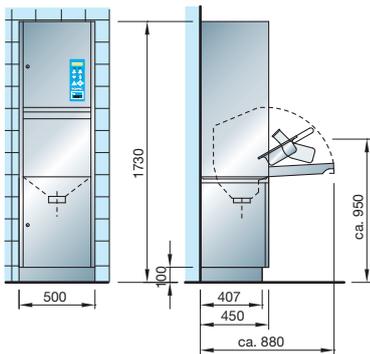

MEIKO TOPIC small combined care unit

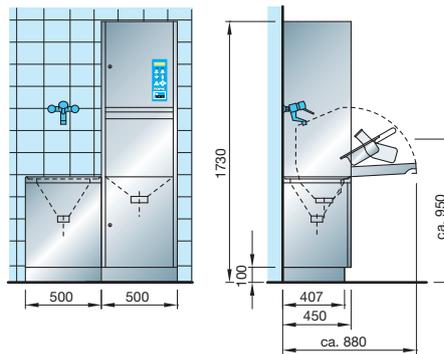

MEIKO TOPIC cabinet model

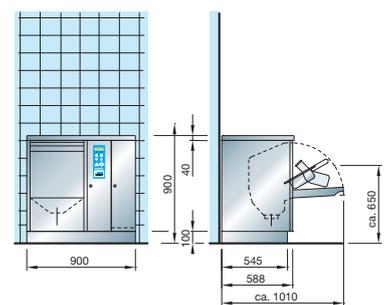

		TOPIC	
		400 V	230 V
Total connected load	kW	4.6	2.8
Feed water connection, cold water		G ¾"	
Feed water connection, hot water 45-60 °C		G ¾"	
Flow pressure	bar	1-5	
Flow	l/min	18	
Drain		DN 100	
Noise level	dB(A)	48.2	



The clean solution



Manufacturer:
MEIKO Maschinenbau GmbH & Co. KG · Englerstr. 3 · 77652 Offenburg, Germany
Phone +49 781 203-0 · info@meiko-global.com · www.meiko.com

760.125.04.06.25/EN

We reserve the right to amend specifications
as part of our product improvement process